



Giving women the freedom to pass urine conveniently



For women with restricted mobility,
continence problems and physical disabilities
as well as the elderly and infirm.

CE Marked
Medically approved
Available on NHS prescription
FDA Compliant

www.whizaway.com

The contemporary solution to an age-old female problem

The Whiz is a small, ingenious device created to revolutionise every woman's life. Hygienic, comfortable and non-invasive, the Whiz means that women can now pass urine in all those situations where previously it has been difficult or embarrassing to do so.

With the Whiz, women with mobility restrictions, disabled or have continence problems, can now pass urine standing up or sitting down, in a wheelchair or car - without undressing. With the Whiz hydration is never a concern even when there are no loos.


Easy to use, anti-bacterial & liquid repellent the Whiz means there is no fuss, no embarrassment and no limits - women are now free to enjoy a normal, active life and stop organising their lives around their bladders.

“This product has made such a difference to my life, I cannot begin to tell you. I have got my confidence back. I am not permanently looking for a loo and it makes me feel good. Wow is what I say. Fantastic.”

Jay T., Birmingham

With a bad back and loos that are too low the Whiz is: “An answer to my prayers”.

Mary F, Oxford

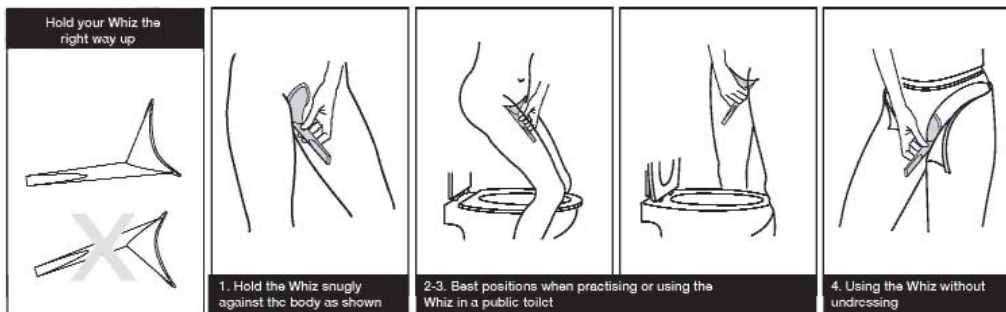

✓ Checklist

- No need to fully undress
- Can be used standing or sitting, indoors or outdoors
- Can be used in restricted spaces such as a car, small toilet etc.
- Small & discreet - easily fits into bag or glove compartment
- Comfortable - made of soft plastic, moulded to the female form
- Designed not to over-fill, spill or splash back
- Self-cleansing and re-usable - so inexpensive & eco-friendly
- Hygienic and hydration friendly
- 'Lily' shape clinically tested by over 1,400 women
- Used by thousands of women including disabled & Mobility impaired
- CE marked, medically approved & anti-bacterial
- Available on NHS Prescription

“With a connecting tube I am able to use the Whiz in the wheelchair and also in the car. It is very quick and empties cleanly into a bottle with no mess. I am very pleased to have solved the problem of how to go to the loo when toilets are inaccessible or locked up. Thank you very much”

Rosemary L, London

Using The Whiz[®]

These are general instructions designed to act as a guide, but remember that every woman is different with individual needs, so adjust these instructions to suit yourself.

The Whiz is a wonderful solution to most urinary problems. However, we are aware that if you have mobility or continence problems you may need to use the Whiz while seated or lying down, in a wheelchair or in a car, in restricted places or difficult situations - all areas where you need to manipulate the flow very carefully. For these circumstances we offer connecting tubes.

To Attach Tubing to the End of the Whiz

1. Put the tubing onto the end of the Whiz, pushing it about 1cm down the nozzle.
2. Give it a small twist and shake to ensure that it is securely attached before using.

Length of Tubing Required

This is a personal choice and depends on how the Whiz is to be used. For some women, 10 cms is enough but for others, the tubing needs to be longer - for example, some women want to run the tubing down the inside of their trouser legs and others like to keep the bottle on the floor in seated situations and need longer tubing to reach the receptacle. Tubing is available in all lengths but if you have any questions please call our Whiz Line: +44 (0) 1865 240572

Basic Instructions

These instructions apply whether seated or lying down

1. Attach tubing securely to end of the Whiz nozzle.
2. Ensure that the Whiz nozzle and tubing are pointing downwards as gravity will promote the flow of urine.
3. Adjust your clothing to allow use (see Tips from Customers).
4. Push underwear aside and place Whiz in position with one hand.
5. Place end of tube into a bottle or attach to a catheter bag and hold with other hand.
6. As you finish, tip the head of the Whiz forward to ensure no dripping or seepage.
7. If you add tubing also remember to raise the Whiz to expel all urine from the tube.

“Wow, wow, wow – what a product! I’ve been waiting for this all my life. Finally an invention that advances feminine hygiene at a stroke”

Shelley G, Birmingham

“No more waiting, squatting, holding-on or hiding! With the Whiz at my side I’m confidently prepared for any emergency.”

KP, Oxford

“I’ve tried the Whiz out and am amazed at how well it works. I have two very young children and always have a problem getting to the toilet when I’m out with them both in tow. Thankfully I now know I will be able to cope with any situation. Thank you for solving a very difficult problem”

Gillian M, Aberdeen

Specific User Instructions

These instructions should be considered in tandem with the basic instructions.

In a Wheelchair

1. First apply the brakes (to avoid accidents such as wheelchair slipping away from under you).
2. Move footplates to the side and place feet firmly on the ground with legs apart.
3. Move forward on the seat so that you will have a downward angle for the Whiz.

In a Car

1. If using the Whiz in the drivers seat, make sure the seat is pushed back to give more room.
2. Move to the front of the seat to allow the Whiz nozzle to point downwards.
3. If you are unable to hold the receptacle make sure it is securely encased

When Seated Indoors

1. Many chair seats tilt backwards, so it may be necessary to raise the body by sitting on a cushion or pillow.
2. It is easier to use the Whiz while seated on a firm rather than soft chair.
3. Move to the front of the seat to allow the Whiz nozzle to point downwards.
4. If you are in bed and sufficiently mobile, move to the edge of the bed and use as directed above.

When in a Lying Position

In bed or lying back, it is important that there is some angle on the nozzle of the Whiz - this may be achieved by sitting or raising the pelvis on a cushion or pillow or finding any comfortable method to give a downward slant to the nozzle of the Whiz and the tubing. Even if the angle is slight, this will help the gravity take the urine away from the body.

Top Tip

Try your Whiz in the shower or bath before venturing out. Once used quietly at home, you will see how easy it is to use and you will then feel confident, whatever the circumstance

Tips from Customers

I love my Whiz but I found that to gain my confidence using it out and about, I needed to practice with it in my shower cubicle. Once I had mastered how to hold it to suit me, I didn't feel at all worried and now I use it all the time.


If you are wearing trousers with a side opening or with an elasticated waist, pull them down slightly

If you have a zip or button opening, they may not go down far enough, so you may need to take the trouser down a little.

If you direct your urine into a bottle, choose one that is comfortable and easy to hold.

I found it easier to open the door of the car and put my feet outside as the steering wheel got in the way.

I found that I could use it lying on my side, with my top leg slightly over my lower leg. This helps with the downward pointing of the Whiz and the flow of urine.

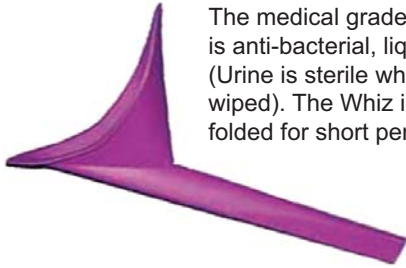

The medical grade Thermoplastic-Elastomer used to make each Whiz is anti-bacterial, liquid repellent, machine washable and easy to dry (Urine is sterile when it leaves the body so it is easily shaken, washed or wiped). The Whiz is also so flexible that it can be temporarily rolled or folded for short periods of time.

The Whiz can be hand-washed, it can also be placed in the dishwasher (or washing machine) in temperatures up to 60°C.

“At my age, I often find myself in difficult circumstances but with the Whiz I have been saved from many embarrassing situations, especially when I am out and about. It has been a real life saver”.

Pat P, Devon

After 2 new hips and 2 new knees following chronic arthritis... “I can go away for the day, tuck it in my bag and wont be stuck for a toilet”.

Jean L, Dundee

“Many thanks – this will probably be the best Crimbo presi ever. I am a wheelchair user and a micro-light pilot. Oh, those small airfield loos are horrible and the flying suit – one piece, plus layers of Damart. Ahhhhhh! Thank you, thank you, thank you.”

Dawn H, Bradford

When, Where, Why...

The Whiz will help when:

- You find it difficult getting to a toilet.
- You find sitting down to urinate a problem.
- You need to pass urine frequently or unexpectedly.
- You have urinary problems requiring a device to direct urine flow.
- You have mobility / physical problems making urination a problem.
- You are wheelchair bound or confined to bed.
- You are fed up with needing people to help you.
- You want to regain your social life, independence, freedom & dignity.
- You are in the military and in action with no place to go.

...and Who?

Disabled Drivers - Elderly women - Women with continence problems & reduced bladder control - Women who have mobility restrictions - Children with disabilities or who need assistance to go to the toilet - Women who have experienced stroke - Women amputees - Women who are wheelchair bound - Women with arthritis or who lack strength in their lower limbs - Women with Multiple Sclerosis - Women with small bowels - Military Personnel - Men with penile amputation or retraction

Whiz Freedom Accessories:



Extension Tube (1):

Sometimes you can get caught in restricted places or difficult situations (kayaking, glider piloting) where you need to manipulate the flow very carefully. Cut to any length in 10cm multiples and easily attached
Can be directed into any receptacle and all known catheter bags
Available with the Whiz on NHS prescription



Hand Sanitizers (2):

Using the Whiz is clean, dry and hygienic, but if you prefer to wash your hands after use, then we offer Purell Instant Hand Sanitizer:
No need for water
Kills 99.9% of germs
Leaves hands soft and smelling fresh



Organza Bag (3):

A smart, discreet and practical way of storing your Whiz:
Homemade
Organza fabric
Available in pink or purple



Explorer Bag (4):

A robust bag for travelling, sporting and outdoor activities:
Machine-made
Strong plastic material
Durable Zip



Whiz Bag (5) :

A general all purpose bag usable anywhere:
Machine-made
Strong plastic material
Durable Zip
Space for any extras
Available in pink purple or black



Drainage Bag (6):

Ideal for temporary or permanent mobility impairment.
2 litre plastic drainage bag
1 metre inlet tube with tap



Gift Sets (7):

The whiz is available with a combination of accessories.
Gift Set A: Whiz freedom, 10cm tube and Organza bag
Gift Set B: Whiz freedom, 10cm tube, Organza bag and 15ml hand sanitizer
Gift Set C: Whiz freedom, 10cm tube, Organza bag, 15ml hand sanitizer and drainage bag

Four easy ways to order your Whiz

- Visit our secure online order form at www.whizaway.com
- Call our sales team on +44 (0) 1865 240572
- Fax us on +44 (0) 1865 246565
- Write to us at : JBOL Ltd., 1 Folly Bridge, Oxford OX1 4LB

Your Whiz on Prescription*

- The Whiz (PIP Code: 304-2827 NHS Order Code: ATD-2002.1)
- The Whiz & Extension Tube (PIP Code: 314-5554 NHS Order Code: ATD-2004.10)

* If you think the Whiz could help your situation, ask your doctor more about the Whiz

Patents Granted: GB2386257, GB2396819, GB23855276
UK Reg. Design: 2100623, 2103229
International Patent Application: 01/74275, PCT/GB2003/000472
Further, patents pending
Manufactured by:
JBOL Ltd, 1 Folly Bridge, Oxford, OX1 4LB, UK
Tel: +44 (0)1860 240572 Fax: +44 (0)1855 246555
(Web: www.whizaway.com / Email: info@whizaway.com)
The Whiz is a registered name & trademark of Jbol Ltd

Distributed by:

JBOL Ltd.
1 Folly Bridge, Oxford, OX1 4LB, UK
Tel: +44 (0) 1865 240572